

ZYKLUSZEIT

00:37:09

H : MIN : SEK

STERILISATIONSDAUER

03:30

MIN : SEK

STERILISATION TEMP.

MIN 134.5 MAX 135.1

°C

F0

118.7

LETZTE TESTS

BOWIE-DICK TEST -


LEAK TEST BESTANDEN
17.03.2026
06:17:41
Nr. 13539

RACK-ID

-

PROGRAMM-NR. / VERSION
5 / 07.08.2020

PROGRAMM-NAME
Bowie-Dick-Test

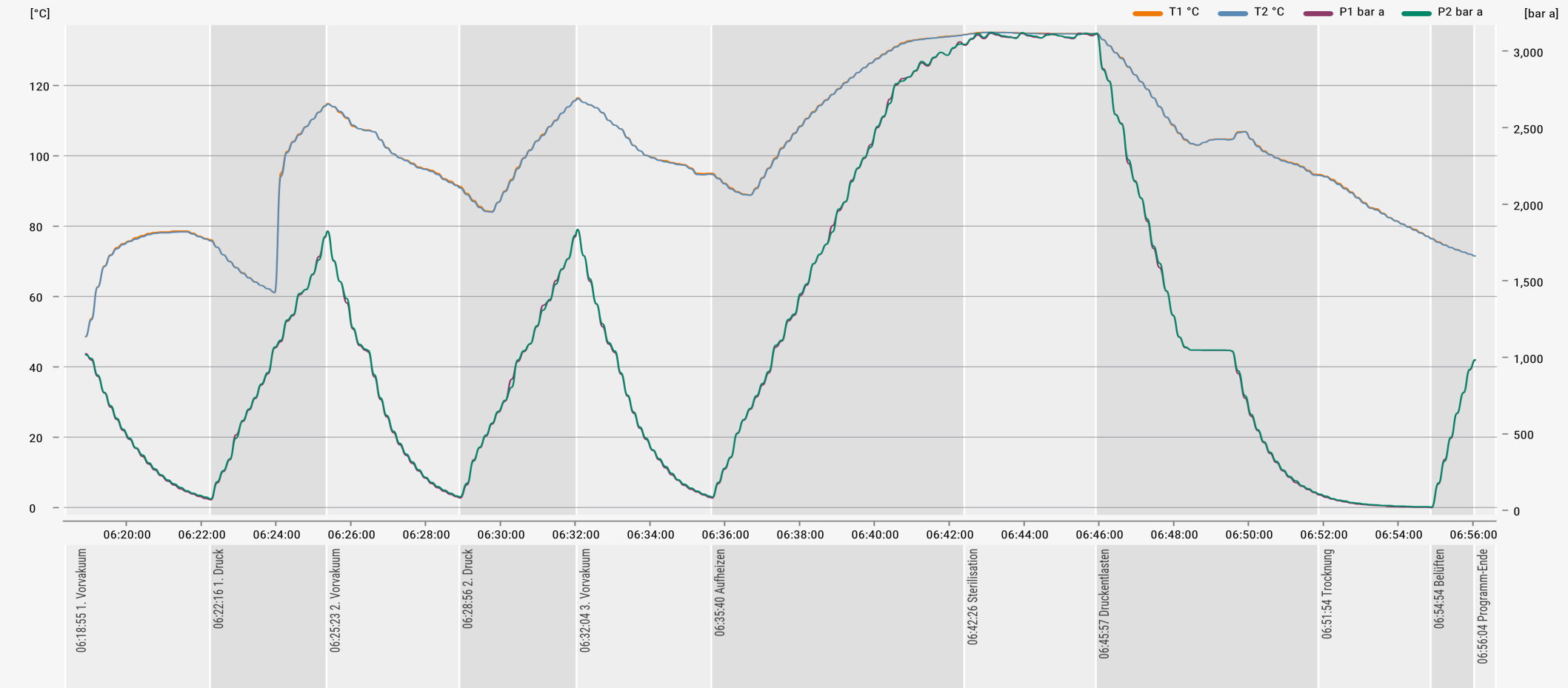

OPERATOR

CHARGE

BEMERKUNGEN ZUR BEDINGTEN FREIGABE

SIGNATURE

User

freigegeben ☐ ja ☐ nein ☐ bedingte Freigabe

Zeit	Meld. Nr.	Meldung	Phase	T1 °C	T2 °C	P1 bar a	P2 bar a	P3 bar a	P4 bar a
06:18:55			1. Vorvakuum	48.50	48.60	1024	1016	3064	4045
06:22:08			1. Vorvakuum	76.50	76.30	78	86	3198	4201
06:22:16			1. Druck	76.10	75.80	70	76	3183	4190
06:25:19			1. Druck	114.30	114.10	1783	1789	3183	4178
06:25:23			2. Vorvakuum	114.90	114.70	1824	1822	3175	4099
06:28:50			2. Vorvakuum	91.70	91.50	89	94	3188	4144
06:28:56			2. Druck	91.30	90.80	80	88	3180	4148
06:32:00			2. Druck	115.70	115.60	1786	1797	3183	4443
06:32:04			3. Vorvakuum	116.50	116.20	1834	1834	3175	4347
06:35:36			3. Vorvakuum	95.00	94.80	84	92	3188	4051
06:35:40			Aufheizen	95.00	94.70	80	86	3184	4365
06:42:17			Aufheizen	134.30	134.20	3060	3046	3206	4158
06:42:26			Sterilisation	134.50	134.40	3039	3044	3181	4201
06:45:49			Sterilisation	134.80	134.70	3101	3109	3080	4292
06:45:57			Druckentlasten	134.80	134.70	3110	3116	3145	4209
06:51:46			Druckentlasten	94.80	94.50	116	121	3188	4342
06:51:54			Trocknung	94.70	94.40	100	106	3179	4407
06:54:47			Trocknung	77.20	77.10	20	24	3195	4216
06:54:54			Belüften	76.50	76.40	19	23	3185	4223
06:55:55			Belüften	72.00	72.00	917	922	3190	4151
06:56:04			Programm-Ende	71.50	71.50	981	984	3181	4147
06:56:04			Programm-Ende	71.50	71.50	981	984	3181	4147